


# WhiteLase SC Series

Blue & UV enhanced Supercontinuum Fiber Lasers

ULTRA BROADBAND  
WHITE-LIGHT LASER

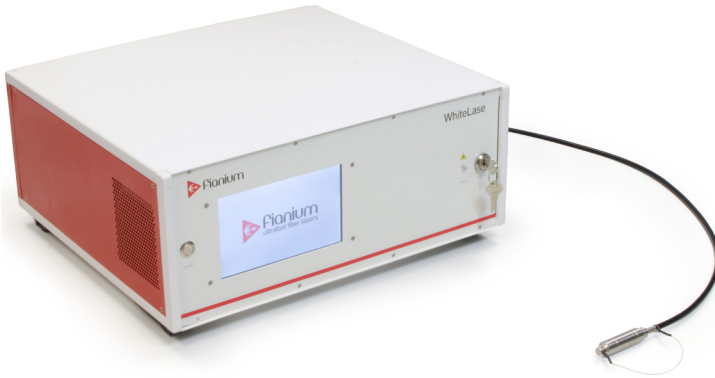



## KEY FEATURES

- Blue or UV enhanced spectrum on all systems
- Fixed or variable repetition rate
- Highest available output power to >20W
- Ultimate visible brightness (350-750nm)
- Touchscreen control with intuitive operation
- Single spatial mode across the output spectrum
- Advanced triggering options
- Modular and upgradable design
- High reliability and low cost of ownership
- Wide range of Plug-and-Play filter options

The **WhiteLase™ SC** range of systems is high-power fibre lasers generating ultrafast broadband supercontinuum light. Operating in the MHz repetition rate range, with picosecond pulses, these systems can be utilised effectively for both steady-state and lifetime measurement. The inherently robust all-fibre design provides unsurpassed performance combined with high reliability and ease-of use. The touchscreen interface enables one-touch access to all laser settings, pre-sets, diagnostics and live system status.

**WhiteLase™ SC-UV;** UV-enhanced system providing a short cut-in wavelength with high brightness from approximately 390nm.

**WhiteLase™ SC400;** Versatile supercontinuum laser with a spectral range from approximately 400nm for visible and near-infrared applications.

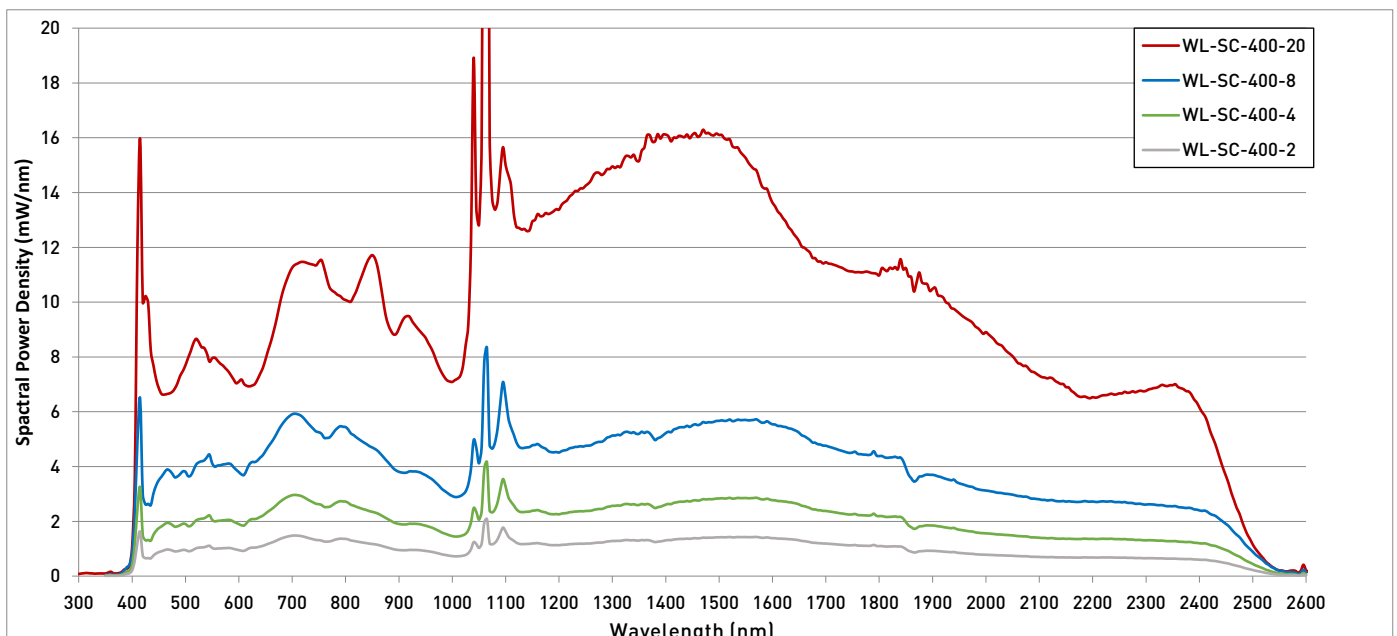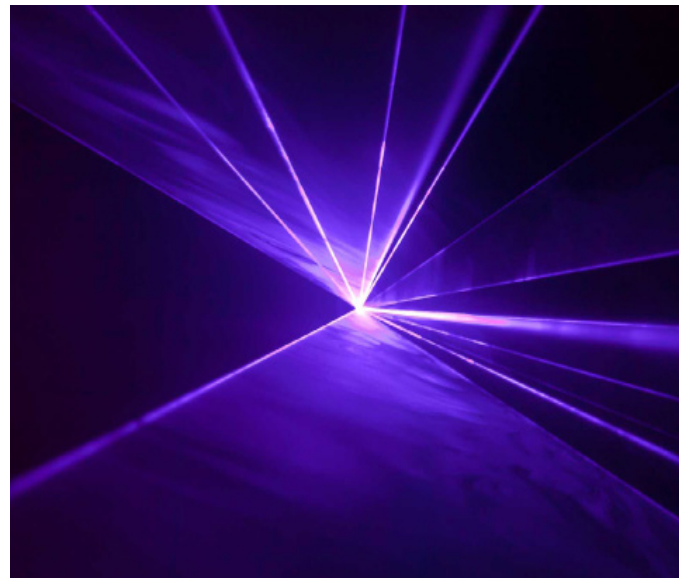

# STANDARD SPECIFICATIONS

	WhiteLase Short Wavelength Supercontinuum Systems				
Model	WL-SC-UV-3	WL-SC-400-2	WL-SC-400-4	WL-SC-400-8	NEW: WL-SC-400-20
Minimum Wavelength	<390nm	<400nm			
Maximum Wavelength	≈2400nm				
Total Power (full spectrum)	>3W	>2W	>4W	>8W	>20W
UV/Visible Power (350-750nm)	>300mW	>300mW	>600mW	>1200mW	>2700mW
Extended Vis Power (350-850nm)	≈450mW	≈450mW	≈900mW	≈1800mW	≈4000mW
Average Spectral Power Density	>1.5mW/nm	>1mW/nm	>2mW/nm	>4mW/nm	>7mW/nm
Fundamental Repetition Rate	40MHz	20MHz	40MHz	80MHz	≥80MHz
Adjustable Repetition Rate Range (with optional pulse picker)	100kHz - 40MHz	100kHz - 20MHz	100kHz - 40MHz	100kHz - 80MHz	N/A
Power Stability	<1%				
Fundamental Pulsewidth	≈6ps				
Output Optic	φ16 x 50mm Collimator				Divergent SM05 compatible
Beam Diameter	≈1.5mm @ 530nm ≈2mm @ 633nm ≈3mm @ 1100nm				N/A
Armoured Fiber length	1.5m				
State of polarisation	Unpolarised				
User Interface	1. Integrated touchscreen graphical user interface 2. PC via USB interface				
Sync (trigger) Outputs	1. NIM Compatible trigger with adjustable delay 2. Oscillator monitor photodiode				
NIM Trigger Specifications	>10ns adjustable delay in 10ps steps <10ps timing jitter				
Cooling	Integrated air cooling				
Dimensions (mm)	19" benchtop chassis - 450 x 390 x 180				

## QUALITY AND RELIABILITY

Fianium have over 10 years' experience designing, high power supercontinuum fiber lasers and have supplied over 1000 units to research and industry. The systems are qualified by industry partners based on:

- 100s of systems tested over 1000s of hours
- Ongoing Highly Accelerated Lifetime Testing (HALT)
- Modules rated to over 20,000 hours Mean Time To Failure (MTTF)

Fianium's unique technology enables leading performance without sacrificing reliability or lifetime. Each system supplied is tested on a module and system level for a combined 500 hours and is backed by an unlimited hours warranty.



All Fianium lasers are manufactured under our Quality Management System in accordance with the ISO 9001:2008 standard.

**FIANIUM UK LTD.**

Tel: +44 2380 458776

Email: [info@fianium.com](mailto:info@fianium.com)

**FIANIUM US INC.**


Tel: +1 541 343 6767

Email: [sales@fianium.com](mailto:sales@fianium.com)

**FIANIUM ASIA LTD.**

Tel: +852 2607 4236

Email: [asia@fianium.com](mailto:asia@fianium.com)

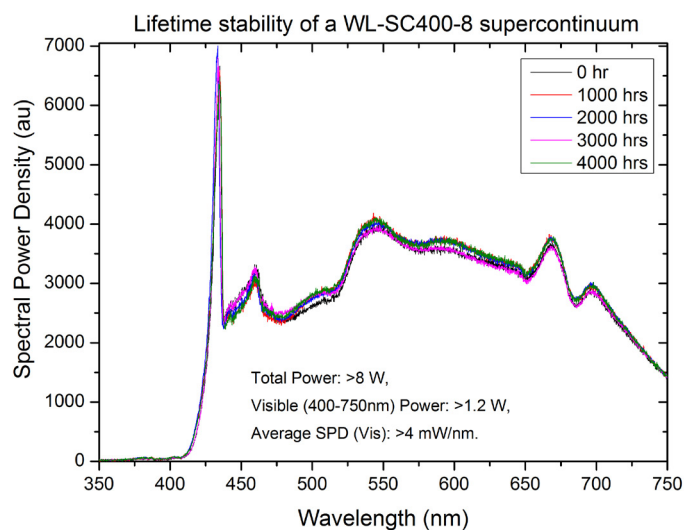

Information contained herein is deemed to be reliable and accurate. Product modification, combination with other products, or use in a specific application may require licensing of 3rd party intellectual property (IP). Customers/users are solely responsible for identifying any such applicable 3rd party IP and obtaining any required licenses or rights. No warranty is made - the customer/user assumes all liability for any infringement of such 3rd party IP. Fianium reserves the right to change the design, specification etc of the products at any time without notice.

